



The Original BioBall® System

Enhanced Flexibility in Primary and
Revision Hip Arthroplasty



BioBall® System
Merete® Innovative Hip Solutions



Merete stands for expert orthopaedic and trauma surgery solutions.

The market for bone surgical medical devices relies upon tested solutions. By developing simple solutions for difficult problems, Merete has blazed trails that have now become gold standards in medical technology. The story of BioBall® as a simple yet unparalleled system has been continued through other Merete products and solutions. Merete GmbH's impressive product families are sophisticated modular systems that allow optimum results under practically any intraoperative circumstances.

A handwritten signature in blue ink, appearing to read 'A. Anapliotis'.

Alexia Anapliotis,
CEO Merete GmbH



Merete GmbH is certified in accordance with EN ISO 13485.
Our Quality Management System is regularly audited by our notified body MedCert (Hamburg).

Table of Contents


4-5	Overview of the BioBall® System
6-11	BioBall® Adapter System
6-7	Overview
8-9	Implant Ordering Information
10-11	Instrument Ordering Information
12-13	BioBall® AdapterSelector™
12	Overview
13	Handling – Step by Step
14-17	BioBall® AdapterSelector™ for Customer Tapers
14-15	Overview
16	Ordering Information
18-22	BioBall® MaxiMotion™ Cup
18-19	Overview
20	Implant Ordering Information
21-22	Instrument Ordering Information
23	BioBall® References

Overview of the BioBall® System

**BioBall®
Metal Head**



**BioBall DELTA™
Ceramic* Head**

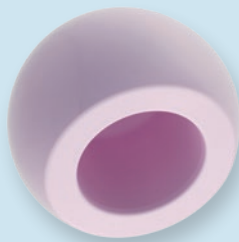

BioBall® Bipolar Duo Head
with preassembled
BioBall® Metal Head



BioBall® AdapterSelector™
Instrument for Intraoperative
Inspection of Taper Geometry




**BioBall® Adapter
Standard 12/14**



**BioBall® Adapter
Standard 14/16**




BioBall® MaxiMotion™ Cup
cemented

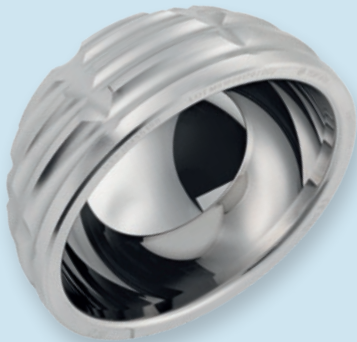


BioBall® MaxiMotion™ Cup
TPS Coating with BONIT®, non-cemented



BioBall® MaxiMotion™ XPE Inlay
with preassembled BioBall® Metal Head

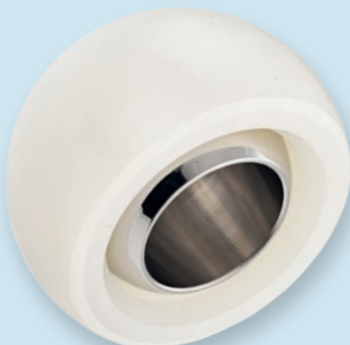

BioBall® MaxiMotion™ XPE Inlay
with preassembled BioBall DELTA™
ceramic* head

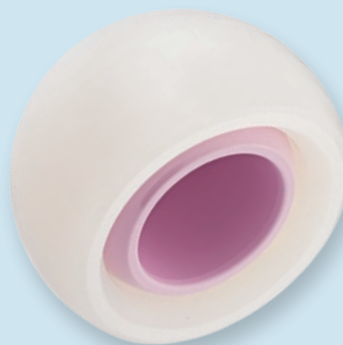

BioBall® Adapter
Offset 12/14



BioBall® Adapter
Offset 14/16



Additional adapter sizes
(tapers / angles) available
on request.


* **Material:**
BIOLOX® delta ceramic from
CeramTec GmbH.



BioBall® Adapter
Standard and Offset



BioBall®
Metal Head

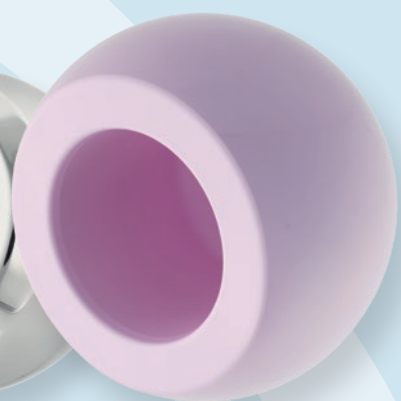

BioBall DELTA™
Ceramic Head

BioBall® Adapter System

The Gold Standard in Revision Surgery

Merete brought the BioBall® System onto the market as a “modular joint prosthesis system” at the end of the 1990s. In the meantime, it has become the gold standard in revision hip arthroplasty surgery. The BioBall® Adapter, consisting of titanium alloy, allows intra-operative correction of neck length as well as antetorsion/retrotorsion and lateralisation/medialisation on in situ stems. This significantly improves the gait pattern and reduces the dislocation risk. With its offset components and special tapers, the system ought to be available in every hospital as a solution for unexpected situations in primary endoprosthetic care as well.

When it comes to revision surgeries, the BioBall® System offers another unique feature: With stems left in situ with previously used tapers, the BioBall® Adapter compensates mild deformation of the taper surface.

Depending on the specific model, BioBall® Adapters are available in sizes ranging from S - 5XL, as standard or offset versions, for 12/14 and 14/16 tapers. Custom adapters for additional tapers are available on request.

Characteristics

- Sliding pair revisions
- Intraoperative correction of neck length
- Intraoperative correction of retro-/antetorsion
- Intraoperative correction of lateralisation/medialisation
- Adjustment of leg length discrepancy in the context of soft tissue management



Scan the QR code and learn more about surgical techniques.

With kind support from
PD Dr. med. Patrick Weber, LMU Munich.

Implant Ordering Information



BioBall® Adapter Standard 12/14 Sterile

Neck length	S (-3.0)	M (0)	L (+3.5)	XL (+7.0)	2XL (+10.5)	3XL (+14.0)	4XL (+17.5)	5XL (+21.0)
Ref.	HM30121	HM30122	HM30123	HM30124	HM30125	HM30126	HM30127	HM30128



BioBall® Adapter Offset 12/14 Sterile

Neck length	M (0)	L (+3.5)	XL (+7.0)	2XL (+10.5)	3XL (+14.0)	4XL (+17.5)	5XL (+21.0)
Offset (mm)	1.1	1.2	1.3	1.5	2.0	2.5	3.0
Ref.	HM30222	HM30223	HM30224	HM30225	HM30226	HM30227	HM30228



BioBall® Adapter Standard 14/16 Sterile

Neck length	M (0)	L (+3.5)	XL (+7.0)	2XL (+10.5)	3XL (+14.0)	4XL (+17.5)	5XL (+21.0)
Ref.	HM30142	HM30143	HM30144	HM30145	HM30146	HM30147	HM30148




BioBall® Adapter Offset 14/16 Sterile

Neck length	2XL (+10.5)	3XL (+14.0)	4XL (+17.5)	5XL (+21.0)
Offset (mm)	1.4	1.5	2.0	2.5
Ref.	HM30445	HM30446	HM30447	HM30448

Implant Ordering Information

BioBall DELTA™ Ceramic Head, suitable for all BioBall® Adapters

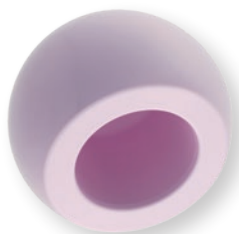

Material:
BIOLOX® delta ceramic*

Size (mm)	Ref.
Ø 28	HM50028
Ø 32	HM50032
Ø 36	HM50036

BioBall® Metal Head suitable for all BioBall® Adapters



Material:
Vivium®**

Size (mm)	Ref.
Ø 28	HM30028
Ø 32	HM30032
Ø 36	HM30036

BioBall® Bipolar Duo Head with preassembled BioBall® Metal Heads, suitable for all BioBall® Adapters



Material:
Vivium®**, UHMWPE



Description	Ref.
Instrument Tray	HM20500


Size (mm)	Ref. Duo Head	Ref. Trial head	Size metal head (mm)
Ø 42	HM30342	HM40342	Ø 28
Ø 43	HM30343	HM40343	Ø 28
Ø 44	HM30344	HM40344	Ø 28
Ø 45	HM30345	HM40345	Ø 28
Ø 46	HM30346	HM40346	Ø 28
Ø 47	HM30347	HM40347	Ø 28
Ø 48	HM30348	HM40348	Ø 28
Ø 49	HM30349	HM40349	Ø 28
Ø 50	HM30350	HM40350	Ø 32
Ø 51	HM30351	HM40351	Ø 32
Ø 52	HM30352	HM40352	Ø 32
Ø 53	HM30353	HM40353	Ø 32
Ø 54	HM30354	HM40354	Ø 32
Ø 55	HM30355	HM40355	Ø 32
Ø 56	HM30356	HM40356	Ø 32
Ø 57	HM30357	HM40357	Ø 32
Ø 58	HM30358	HM40358	Ø 32

*BIOLOX® delta is a registered trademark of CeramTec GmbH.

**Vivium® is a registered trademark of Merete GmbH.

Instrument Ordering Information

Instrument Tray

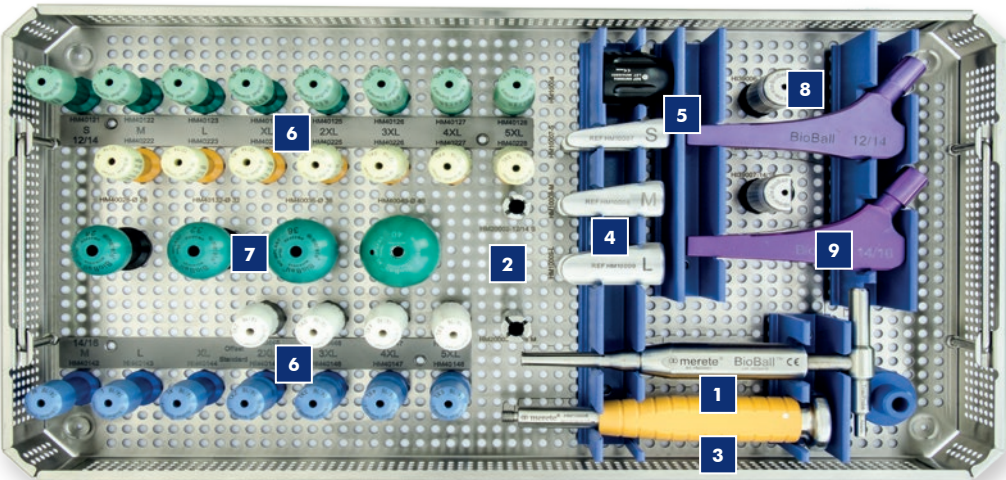

Description	Ref.
Instrument Tray	HM30770

1 Separator

Description	Ref.
Separator	HM20001



Adapter Sleeve for Adapter 12/14 S	Ref.
	HM20002
Adapter Sleeve for Adapter 14/16 M	Ref.
	HM20003

3 Universal Handle

Description	Ref.
Universal Handle	HM10005



Separator Wedge	Size	Ref.
	S	HM10007
	M	HM10008
	L	HM10009
Head Impactor		Ref.
		HM10004

Instrument Ordering Information

6 Trial Adapter



Length	Ref. Standard 12/14	Ref. Offset 12/14	Ref. Standard 14/16	Ref. Offset 14/16
S (-3.0)	HM40121	—	—	—
M (0)	HM40122	HM40222	HM40142	—
L (+3.5)	HM40123	HM40223	HM40143	—
XL (+7.0)	HM40124	HM40224	HM40144	—
2XL (+10.5)	HM40125	HM40225	HM40145	HM40445
3XL (+14.0)	HM40126	HM40226	HM40146	HM40446
4XL (+17.5)	HM40127	HM40227	HM40147	HM40447
5XL (+21.0)	HM40128	HM40228	HM40148	HM40448


7 BioBall® Trial Heads



Ø	Ref.
28	HM40028
32	HM40132
36	HM40036
40	HM40040

Additional sizes on request.

X-ray Template

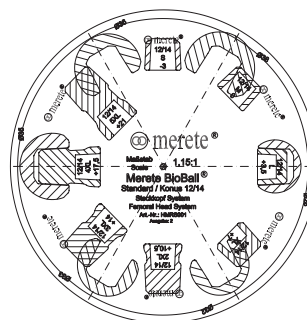

For BioBall® Adapter	Ref.
12 / 14 Standard	HMRS0001
12 / 14 Offset	HMRS0005
14 / 16 Standard	HMRS0002
14 / 16 Offset	HMRS0006

8 BioBall® AdapterSelector™



For taper	Ref.
12/14	HI39006
14/16	HI39007

9 Offset PositionAssistant



Name	Ref.
Offset PositionAssistant 12/14	HM39106
Offset PositionAssistant 14/16	HM39107

BioBall® AdapterSelector™

How do you identify and why do you document taper geometry on in situ stems during revision?
Four reasons why you should know the BioBall® AdapterSelector™.

- 1 Once a surgeon has decided to leave an existing prosthesis stem in place during revision, visual and haptic inspection of the smooth, reflective surfaces is often not enough to determine what the taper is made of. As a technical and mechanical testing instrument, the BioBall® AdapterSelector™ provides information on whether this is the previously defined taper and whether it is damaged.
- 2 Many manufacturers offer hip stems with different taper geometries. In addition, patients from other countries, or those who underwent surgery abroad, often have very old or unfamiliar models and no endoprosthesis record cards. The patented BioBall® AdapterSelector™ helps surgeons inspect the stem taper to determine the correct BioBall® Adapter with great certainty.
- 3 Documented proof that an intraoperative fit check was performed also offers additional security from a legal perspective. If you do a check using the AdapterSelector™ prior to using the BioBall® System, you can document that check in your surgical report.
- 4 The BioBall® AdapterSelector™ is the only testing instrument worldwide approved for testing taper geometry (incl. FDA). No other instrument in the world allows you to perform an approved, recognised taper geometry check and thus ensure that your selected BioBall® Adapter will fit properly.



Scan the QR code
and learn more
about handling the
BioBall® AdapterSelector™.



Handling - Step by Step

Step 1

Remove the existing head of the in situ stem.



Step 2

The taper should be clean and dry before the BioBall® AdapterSelector™ is inserted.



Step 3

Apply slight pressure and turn clockwise to place the BioBall® AdapterSelector™ onto the stem taper. Check whether the taper's flat face is positioned between the two arrows.

If it is positioned above or below these markings on the BioBall® AdapterSelector™, the stem taper is not the same as the taper indicated on the BioBall® AdapterSelector™.



Step 4

Then visually inspect the lateral accuracy of fit and check to see whether there is a gap in the upper or lower taper region between the stem taper and the BioBall® AdapterSelector™.



Step 5

Test the clamp connection on the BioBall® AdapterSelector™ with a sideways tipping movement. If this movement causes rattling or tipping in the BioBall® AdapterSelector™, the stem taper is not the same as the taper indicated on the BioBall® AdapterSelector™.



Step 6

After checking the taper geometry with the help of the BioBall® AdapterSelector™, examine the entire taper surface. Start with the taper's flat face visible in the opening of the BioBall® AdapterSelector™. Then remove the BioBall® AdapterSelector™ and examine the entire taper surface.



BioBall® AdapterSelector™ for Custom Tapers – the Flexible-Use Head System for Various Taper Geometries.

There are still no uniform standards in place for prosthesis stem tapers. Implant manufacturers continue to use tapers designed to their own individual specifications that differ from the others in terms of geometry structure and surface. Neck length sizes (S, M, L, etc.) are not standardised and may vary significantly from manufacturer to manufacturer.

Besides standard and offset **BioBall® Adapters**, we also offer **BioBall® Adapters for custom tapers**, for a wide variety of taper geometries from other manufacturers.

Characteristics

BioBall® Adapter for Custom Tapers

The adapter can compensate misalignment between the prosthesis stem and the acetabular cup. Neck length and offset can be individually selected and adjusted.

- Sliding pair revisions
- Intraoperative correction of retro-/antetorsion
- Intraoperative correction of neck length
- Intraoperative correction of lateralisation/medialisation
- Adjustment of leg length discrepancy in the context of soft tissue management



**BioBall® AdapterSelector™ MST1
and BioBall® AdapterSelector™ MSV4**

For final verification of the taper of the existing stem, and to ensure conformity with legal documentation requirements, always use the recommended BioBall® AdapterSelector™. This is the only reliable means of selecting the appropriate BioBall® Adapter and guaranteeing its technical and medical functionality.

Selection of different Taper Variants

Manufacturer	Name	Taper	Checking Tapers with the BioBall® AdapterSelector™							
			12/14	14/16	MST1	MSZI	MSSR	MSBG	MSV4	MSPC
Biomet / Zimmer ²	12/14	12/14	x							
	Type I	11/13			x					
	6 Degree Taper	10/12				x				
DePuy ³	Articul/eze® Taper ²	12/14	x							
	Large Taper	14/16		x						
	S-ROM® Taper ²	11/13					x			
Waldemar Link ⁴	12/14	12/14	x							
	14/16	14/16		x						
Smith & Nephew ⁵	12/14	12/14	x							
Stryker ⁶ / Howmedica ⁷	C-Taper	12/14	x							
	14/16*	14/16		x						
	V40™	11/12							x	
	6° Taper	14/16						x		
	PCA® Taper ⁶	13/14								x

BioBall® Adapters for custom tapers are only approved for use in combination with BioBall® Metal Heads. No biomechanical testing information is available on the usage of BioBall® Adapters with hip stems from other manufacturers. Consequently, only manufacturer-approved extensions may be used.

* Please see the overview for information on which BioBall® Adapters are available for custom tapers. Our experienced staff would be glad to provide you with more in-depth consultation.



² The name **Zimmer** is a registered trademark of Zimmer, Inc., Warsaw Ind., US/Biomet is a registered trademark of BIOMET Inc., Warsaw Ind., US

³ The names **DePuy** and **S-ROM** are registered trademarks of DePuy Synthes, Inc. Warsaw Ind., US

⁴ The name **Waldemar Link** is a registered trademark of Waldemar Link GmbH & Co. KG, 22339 Hamburg, Germany

⁵ The name **Smith&Nephew** is a registered trademark of Smith&Nephew Plc, WC2N 6LA, London, GB

⁶ The name **Stryker** is a registered trademark of Stryker Corp., Kalamazoo, MI US

⁷ The names **Howmedica** and **ABG** are registered trademarks of Howmedica Osteonics Corp., Mahwah, NJ US

Ordering Information

BioBall® AdapterSelector™

Name	Ref.
BioBall® AdapterSelector™ 12/14	HI39006
BioBall® AdapterSelector™ 14/16	HI39007
BioBall® AdapterSelector™ MST1	HI39001
BioBall® AdapterSelector™ MSV4	HI39002

Name	Ref.
BioBall® AdapterSelector™ MSZI	HI39004
BioBall® AdapterSelector™ MSBG	HI39008
BioBall® AdapterSelector™ MSPC	HI39009
BioBall® AdapterSelector™ MSSR	HI39010

BioBall® Adapter Standard for Custom Tapers

BioBall® Adapter	MSZI	MST1	MSV4	MSBG	MSPC	MSSR
Neck length	(10/12)	(11/13)	(11/12)	(14/16)	(13/14)	(11/13)
S (-3mm)	HM33121	—	—	—	—	—
M (0mm)	HM33122	HM36002	HM34122	HM31142	HM31132	HM31152
L (3.5mm)	HM33123	HM36003	HM34123	HM31143	HM31133	HM31153
XL (7mm)	HM33124	HM36004	HM34124	HM31144	—	HM31154
2XL (10.5mm)	HM33125	HM36005	HM34125	HM31145	—	—
3XL (14mm)	HM33126	HM36006	HM34126	—	—	—





BioBall® Adapter Offset for Custom Tapers


BioBall® Adapter	MSZI	MST1	MSV4	MSBG	MSPC	MSSR
Neck length	(10/12)	(11/13)	(11/12)	(14/16)	(13/14)	(11/13)
M (0mm)	—	HM36022	HM34222	—	—	—
L (3.5mm)	—	HM36023	HM34223	—	—	—
XL (7mm)	—	HM36024	HM34224	—	—	—
2XL (10.5mm)	—	HM36025	HM34225	—	—	—
3XL (14mm)	—	HM36026	HM34226	—	—	—

Custom designs are available for other tapers on request.

Extra Trays (as an add-on or as needed; with trial adapters; non-sterile)

Name	Ref.
BioBall® Instrument Tray MST1	HM30730
BioBall® Instrument Tray MSV4	HM30750
BioBall® Instrument Tray MSBG/MSPC/ MSSR	HM30740
BioBall® Instrument Tray with General Instruments	HM30785
BioBall® Instrument Tray MSZI	on request

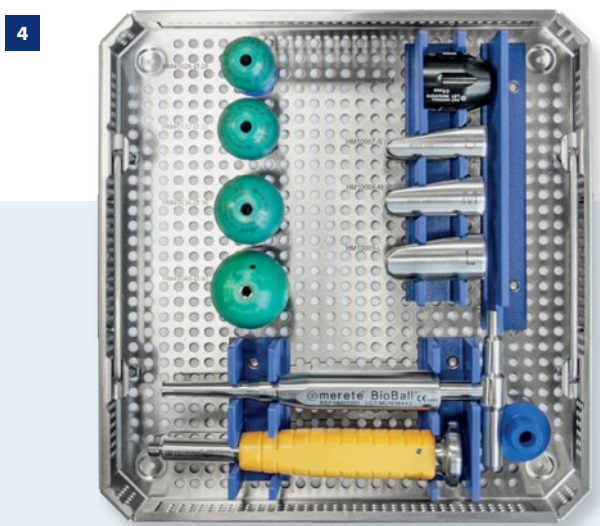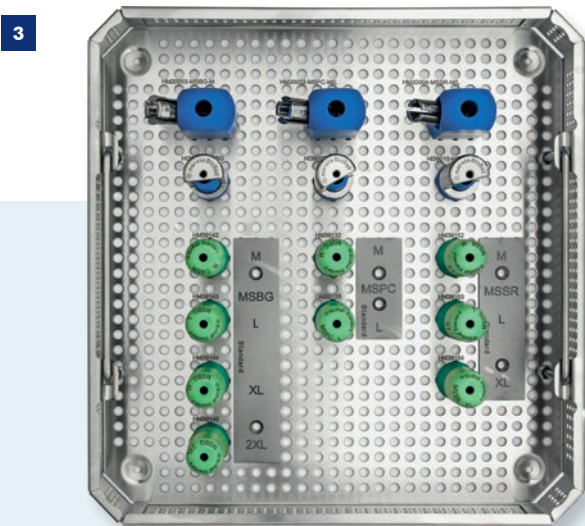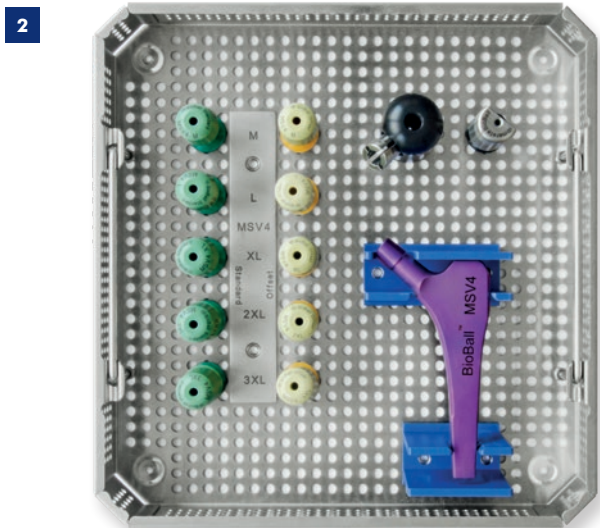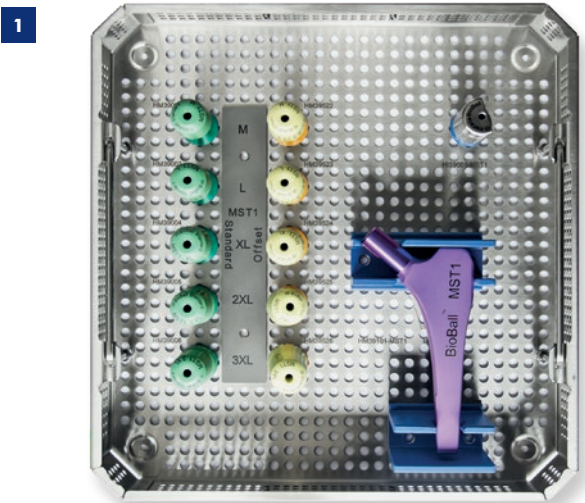

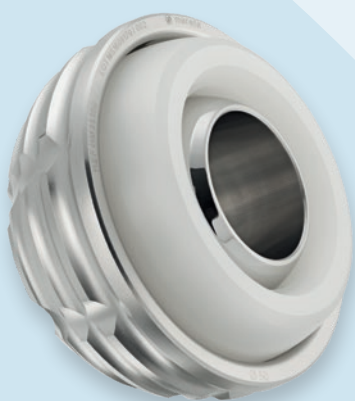

BioBall® MaxiMotion™ Cup,
cemented



BioBall® MaxiMotion™ Cup
TPS Coating with BONIT®,
non-cemented



BioBall® Adapter
Standard

BioBall® MaxiMotion™ Cup

Modular Dual Mobility Cup

The proven BioBall® System can now be combined with the Dual Mobility Concept.


Dislocation is one of the most common complications following primary or revision arthroplasty. Dual Mobility implants improve prosthetic stability and increase range of motion dramatically.

The BioBall® MaxiMotion™ Dual Mobility Cup is a logical expansion to the Merete® BioBall® Adapter System (to be used only in combination with that system).

The result? A great deal of intra-operative flexibility (neck length/offset) plus minimal dislocation risk to give your patients maximum mobility support.

Why the BioBall® MaxiMotion™ Cup?

- The Dual Mobility Concept combines the low friction of small-diameter heads and the exceptional joint stability of large-head prosthetics.
- Effectively lowers dislocation risk by combining the Dual Mobility Concept with BioBall® Adapters.
- Combines with all adapter sizes (standard and offset) in the tried-and-true BioBall® System for even more intraoperative flexibility.
- Available as a cemented or cementless system, with a preassembled 28 mm-diameter BioBall® Metal or BioBall DELTA™ Ceramic Head.
- Wear behaviour tests in an accredited testing laboratory confirm low wear rates.

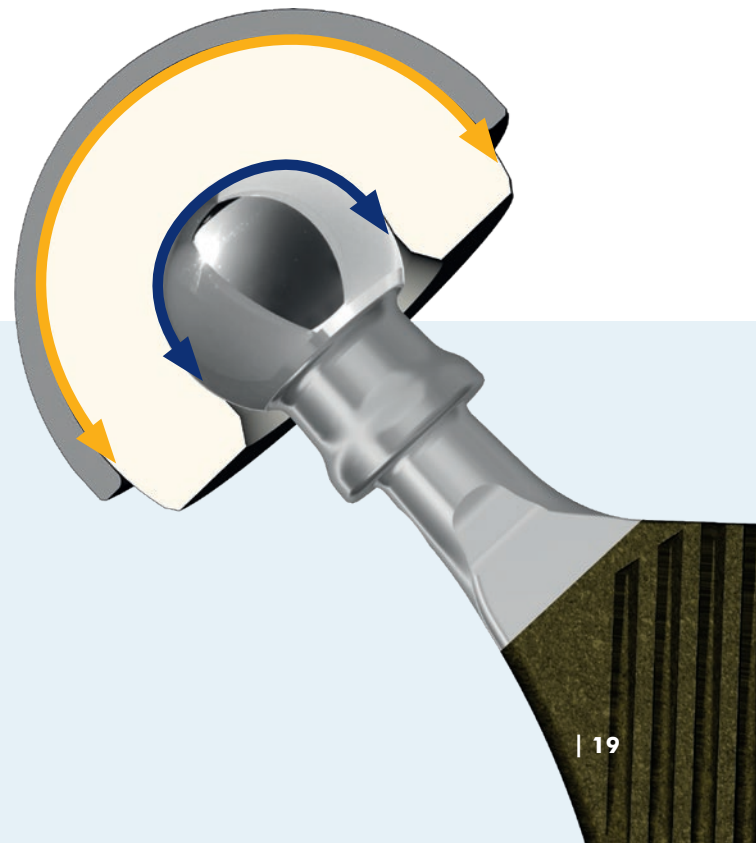

Scan the QR code and learn more about surgical techniques.

With kind support from Head Physician
Christoph Kruis, Rotkreuzklinik Lindenberg.

Implant Ordering Information

BioBall® MaxiMotion™ Cup non-cemented

Material: Vivium™, TPS coating with BONIT®



Size (mm)	Ref.
Ø 46	HM35346
Ø 48	HM35348
Ø 50	HM35350
Ø 52	HM35352
Ø 54	HM35354
Ø 56	HM35356

Size (mm)	Ref.
Ø 58	HM35358
Ø 60	HM35360
Ø 62	HM35362
Ø 64	HM35364
Ø 66	HM35366
Ø 68	HM35368

BioBall® MaxiMotion™ Cup Cemented

Material: Vivium™



Size (mm)	Ref.
Ø 46	HM35146
Ø 48	HM35148
Ø 50	HM35150
Ø 52	HM35152
Ø 54	HM35154
Ø 56	HM35156

Size (mm)	Ref.
Ø 58	HM35158
Ø 60	HM35160
Ø 62	HM35162
Ø 64	HM35164
Ø 66	HM35166
Ø 68	HM35168

BioBall® MaxiMotion™ XPE Inlay*

Material: UHMWPE with BioBall DELTA™ Ceramic Head (28 mm diameter)



Size (mm)	Ref.
Ø 46	HM35669
Ø 48	HM35670
Ø 50	HM35671
Ø 52	HM35672
Ø 54	HM35673
Ø 56	HM35674

Size (mm)	Ref.
Ø 58	HM35675
Ø 60	HM35676
Ø 62	HM35677
Ø 64	HM35678
Ø 66	HM35679
Ø 68	HM35680

BioBall® MaxiMotion™ XPE Inlay

Material: Vivium™, UHMWPE with BioBall® Metal Head (28 mm diameter)



Size (mm)	Ref.
Ø 46	HM35069
Ø 48	HM35070
Ø 50	HM35071
Ø 52	HM35072
Ø 54	HM35073
Ø 56	HM35074

Size (mm)	Ref.
Ø 58	HM35075
Ø 60	HM35076
Ø 62	HM35077
Ø 64	HM35078
Ø 66	HM35079
Ø 68	HM35080

*BIOLOX® delta is a registered trademark of CeramTec GmbH.

Instrument Ordering Information

Instructions for digital planning

Merete® hip products are included in several databases of digital surgical planning tools. For more information about which systems are supported, contact Merete GmbH.


mediCAD®
The Orthopedic Solution

X-ray template

Name	Ref.
For BioBall® MaxiMotion™ non-cemented	HMRS112
For BioBall® MaxiMotion™ cemented	HMRS114



Instrument Tray

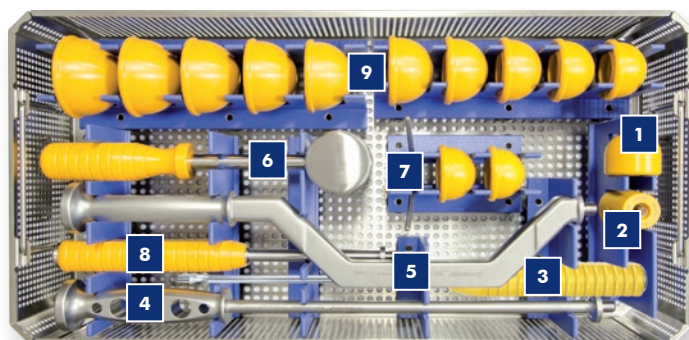

Name	Ref.
Instrument Tray	HM35506

Name	Ref.
1 Impactor Head	HM35508
2 Impactor Tip	HM35505
3 Corrective Impactor	HM35509
4 Seating Instrument (straight)	HM35500
5 Seating Instrument (curved)	HM35501
6 Slotted Hammer	AI00048
7 Aiming Device	HM35502
8 Head Impactor Handle	HI70038

9 Impactor Plate




Size (mm)	Ref.
Ø 46	HM35546
Ø 48	HM35548
Ø 50	HM35550
Ø 52	HM35552
Ø 54	HM35554
Ø 56	HM35556

Size (mm)	Ref.
Ø 58	HM35558
Ø 60	HM35560
Ø 62	HM35562
Ø 64	HM35564
Ø 66	HM35566
Ø 68	HM35568

Reamer Ordering Information

Spherical Reamer Tray

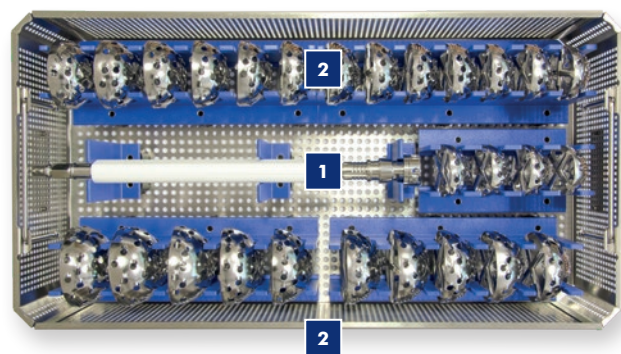


Description	Ref.
Spherical Reamer Tray	HF12080

Description	Ref.
1 Handle for Reamer with AO coupling	HF13006

Optional

Description	Ref.
Handle for Offset Reamer with AO coupling	HF13010

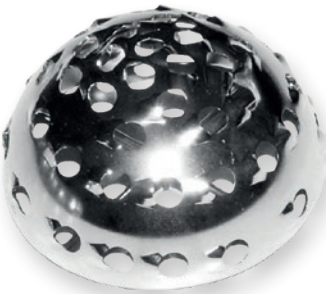

2 Spherical Reamer, individual

Size (mm)	Ref.
Ø 44	HF12044
Ø 45	HF12045
Ø 46	HF12046
Ø 47	HF12047
Ø 48	HF12048
Ø 49	HF12049
Ø 50	HF12050
Ø 51	HF12051
Ø 52	HF12052
Ø 53	HF12053
Ø 54	HF12054
Ø 55	HF12055
Ø 56	HF12056
Ø 57	HF12057

Size (mm)	Ref.
Ø 58	HF12058
Ø 59	HF12059
Ø 60	HF12060
Ø 61	HF12061
Ø 62	HF12062
Ø 63	HF12063
Ø 64	HF12064
Ø 65	HF12065
Ø 66	HF12066
Ø 67	HF12067
Ø 68	HF12068
Ø 69	HF12069
Ø 70	HF12070

BioBall® References

- P. Weber, A. Steinbrück, A. C. Paulus, M. Woiczinski, F. Schmidutz, A. Fottner, V. Jansson (2017). „Gelenkteilwechsel in der Hüftarthroplastik. Was dürfen wir kombinieren?“ Orthopäde DOI 10.1007/s00132-016-3380-4
- Kock, H.J., Cho, C., Buhl, K., Hillmeier, J., Huber, F.X. „Long-term outcome of a modular head adapter system in revision hip arthroplasty for multimorbid patients“ Abstract Nr. 52, Vortrag EHS 2016 - München
- Bloch, B. and S. West (2015). Early results of the BioBall Taper Adaptor in revision total hip arthroplasty. British Hip Society Annual Meeting. London.
- Kock, H. J., C. Cho, F.-X. Huber and J. Hillmeier (2015). 10-year-results After Treatment Of Dislocating Total Hip Arthroplasty Using A Novel Head Adapter System. AAOS.
- Hoberg, M., C. Konrads, S. Huber, S. Reppenhagen, M. Walcher, A. Steinert, T. Barthel and M. Rudert (2015). „Outcome of a modular head-neck adapter system in revision hip arthroplasty.“ Arch Orthop Trauma Surg.
- Friedrich, M. J., S. Gravius, J. Schmolders, M. D. Wimmer and D. C. Wirtz (2014). „Biologisch azetabuläre Defektrekonstruktion beim Hüftendoprothesenwechsel mittels „Impaction Grafting“ und azetabulärem Rekonstruktionsring.“ Operative Orthopädie und Traumatologie 26(2): 126-140.
- Weber, P. and V. Jansson (2014). Teilwechsel an der Hüft. Was dürfen wir kombinieren? Orthopädische Nachrichten. Köln, Biermann Verlag GmbH. 01.2014.
- Helwig, P., L. Konstantinidis, A. Hirschmüller, A. Bernstein, O. Hauschild, N. P. Südkamp and B. G. Ochs (2013). „Modular sleeves with ceramic heads in isolated acetabular cup revision in younger patients-laboratory and experimental analysis of suitability and clinical outcomes.“ Int Orthop 37(1): 15-19.
- Jack, C. M., D. O. Molloy, W. L. Walter, B. A. Zicat and W. K. Walter (2013). „The use of ceramic-on-ceramic bearings in isolated revision of the acetabular component.“ Bone Joint J 95-B(3): 333-338.
- Vaishya, R., M. Sharma and R. R. Chaudhary (2013). „Bioball universal modular neck adapter as a salvage for failed revision total hip arthroplasty.“ Indian Journal of Orthopaedics 47(5): 519.
- Woelfle, J. V., C. R. Fritzl, H. Reichel and D. Wernerus (2013). „Significantly Reduced Leg Length Discrepancy and Increased Femoral Offset by Application of a Head-Neck Adapter in Revision Total Hip Arthroplasty.“ J Arthroplasty.
- Perka, C., B. Fink, M. Millrose, U. Sentürk, M. Wagner, J. Schröder, H. Bail, R. Ascherl, A. Pruss, K. Thiele and C. Götz (2012). Revisionsendoprothetik. AE-Manual der Endoprothetik. L. Claes, P. Kirschner, C. Perka and M. Rudert, Springer Berlin Heidelberg: 441-587.
- Croce, A., M. Ometti and P. Dworschak (2011). „A580. MINIMAL INVASIVE REVISION SURGERY WITH MODULAR NECK ADAPTORS (BIOBALL).“ Journal of Bone & Joint Surgery, British Volume 93-B(SUPP IV): 422-422.
- Kretzer, J. P., R. Sonntag, J. Reinders, E. Jakubowitz, M. Thomsen and C. Heisel (2010). Fretting and Metal Release of Modular Neck Total Hip Arthroplasty. 56th Annual Meeting of the Orthopaedic Research Society. New Orleans.
- Dürr, H. R. (2009). „[The problem of fractures of ceramic heads. What should be done?].“ Orthopäde 38(8): 698-703.
- Kretzer, J. P., E. Jakubowitz, M. Krachler, M. Thomsen and C. Heisel (2009). „Metal release and corrosion effects of modular neck total hip arthroplasty.“ International Orthopaedics 33(6): 1531-1536.
- Fink, B. and A. Grossmann (2008). „[Technique of implantation of a cementless press-fit cup in revisions with severe bone defects].“ Oper Orthop Traumatol 20(2): 157-167.
- Kleffner, B., M. Morlock and R. Schröder (2008). Werkstoff-und designrelevante Aspekte von Pfannenrevisionsimplantaten. Revisionsendoprothetik der Hüftpfanne. D. Wirtz, C. Rader and H. Reichel, Springer Berlin Heidelberg: 201-229.
- Kircher, J., P. Bergschmidt, R. Bader, D. Kluess, E. Besser-Mahuzir, A. Leder and W. Mittelmeier (2007). „Die Bedeutung der Gleitpaarung beim jüngeren Endoprothesenpatienten.“ Der Orthopäde 36(4): 337-346.
- Bader, R., R. Barbano and W. Mittelmeier (2005). „Treatment of recurrent dislocation associated with impingement after revision total hip arthroplasty.“ Acta Orthop Belg 71(1): 98-101.
- Gradingner, R., R. Burkart and M. Goebel (2005). MML-System (ESKA-Implants). Modulare Revisionsendoprothetik des Hüftgelenks. P. Thümler, R. Forst and G. Zeiler. Heidelberg, Springer Medizin Verlag: 258-263.
- Kock, H. J., C. Niewöhner, J. Hillmeier and P. J. Meeder (2004). Frühergebnisse der Behandlung von Hüftprothesenluxationen mit einem modularen Steckkopfsystem bei multimorbiden Patienten. 171. Jahrestagung der Vereinigung Niederrheinisch-Westfälischer Chirurgen. Köln.

Merete GmbH

Alt-Lankwitz 102
12247 Berlin

Tel. +49 (0)30 77 99 80-0
Fax +49 (0)30 76 68 03 61

export@merete-medical.com
www.merete-medical.com