


# Noviomagus Mini Mill



Today's technology guarantees high quality chips

# High quality construction

The Noviomagus Mini Mill is unique because of its superior design. The Mini Mill set only comprises a few parts: a housing, a push block, a screw nut, a handle and various milling drums. All parts of the Noviomagus Mini Mill are machined from stainless steel, which makes it an extremely solid and durable device. There is no risk for abrasion. Additionally, the milling drums are made of a special hardened stainless steel alloy. Even in case of frequent use they do not show signs of wear. This combination guarantees a high quality graft without metal contamination.

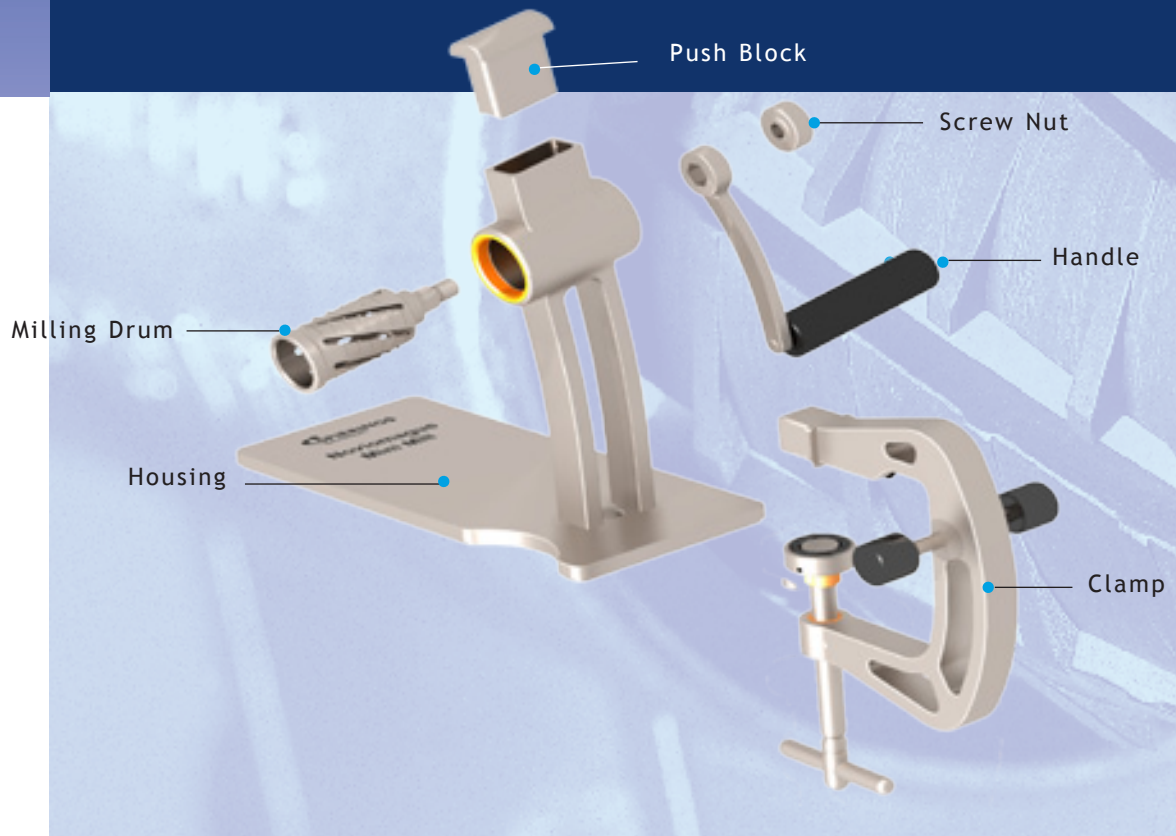

# Easy and optimal use

The teeth of the milling drums feature a unique patented technology. They are spirally positioned and razor sharp. With the stable base plate and optional table clamp, the Noviomagus Mini Mill is easy to handle and it takes little strength to mill bone parts. Because of the optimal design of the housing, the bone is milled without any bone parts left inside the funnel. It can mill bone parts in less than 10 seconds.


*spirally positioned teeth*



# Can be used instantaneously

The Noviomagus Mini Mill guarantees bone chips with an intact trabecular bone structure, because it breaks out chips. The Noviomagus Mini Mill can be used for milling bone harvested from various sites like the iliac crest, jaw and spine. The chips can be used instantaneously during surgery.

The *Extra Fine* milling drum is suitable for spinal and trauma surgery. The *Extra Extra Fine* milling drum is designed to produce extra small bone chips which can be used for maxillo-facial and hand-foot surgery. The *Extra Extra Extra Fine* milling drum produces mini chips for dental surgery.

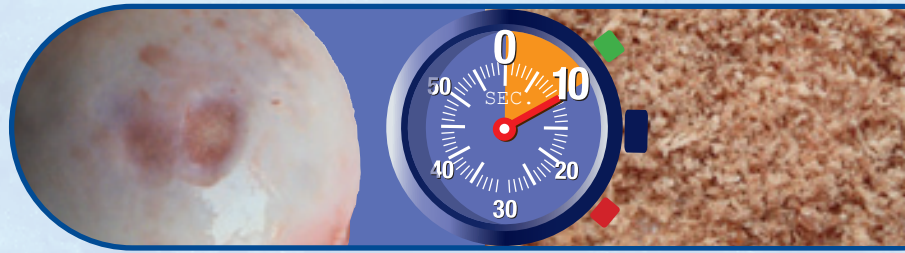

*Less than 10 seconds!*



## Advanced benefits

- Milling of bone particles takes less than 10 seconds
- 3 years guarantee for hospital use
- Three types of milling drums for small chip sizes
  - Extra Fine (3 - 5 mm)
  - Extra Extra Fine (1 - 3 mm)
  - Extra Extra Extra Fine (<1 mm)
- No risk of metal abrasion
- Long life-span without wear of milling drums
- Controlled bone particle size
- Easy cleaning
- Detachable table clamp for additional stability

Today's technology guarantees high quality chips


# Ordering information

## Noviomagus

### Mini Mill

Noviomagus Mini Mill:

Mini Housing	SBM-1002-29
Mini Push Block	SBM-1002-30
Mini Handle	SBM-1002-31
Mini Screw Nut	SBM-1002-32
Noviomagus Mini Milling Drum Extra Fine	SBM-1002-34
Noviomagus Mini Milling Drum Extra Extra Fine	SBM-1002-35
Noviomagus Mini Milling Drum Extra Extra Extra Fine	SBM-1002-33
Sterilisation Tray Mini Mill	SBM-1002-36
Clamp Noviomagus Mini Mill	SBM-1002-37




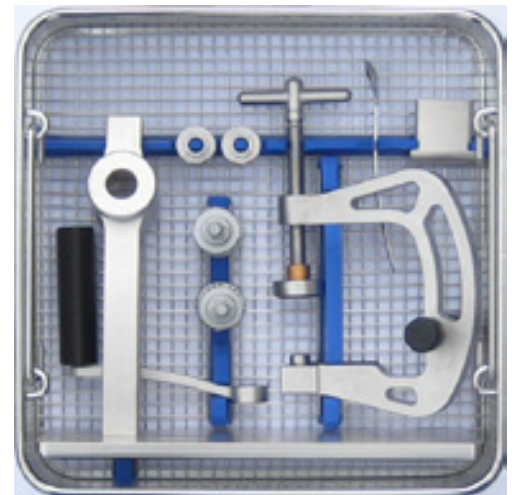


The Noviomagus Mini Mill is CE marked and FDA approved. US Patent No. 6,318,651.

Other Noviomagus products manufactured by Spierings Orthopaedics B.V. are: Noviomagus Impactor Set, Noviomagus Femoral Head Reamer Set, Noviomagus Bone Mill and Noviomagus Revision Meshes.

No rights can be derived from the information in this brochure.

Copyright: Spierings Orthopaedics B.V.

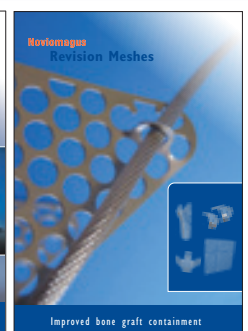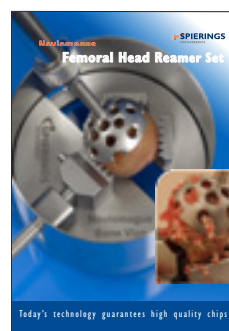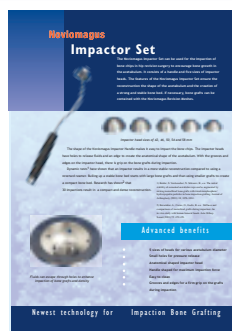

HB-Medical  
 Tivolilaan 205  
 6824 BV Arnhem  
 The Netherlands  
 T: +31 -(0)24-3618804  
 F: +31 -(0)26-3882550  
 E: info@hb-medical.nl  
 I: www.hb-medical.nl