

Noviomagus Micro Mill


The Noviomagus Micro Mill is a compact bone mill designed for implantology procedures, such as a sinus lift. It comprises only a few parts: Micro Housing, Micro Push Block, Micro Handle, Micro Chip Cup, Micro Wrench and two Micro Milling Drums. The milling drums are made of a specially hardened stainless steel alloy, which minimizes wear even in case of frequent use.

The Noviomagus Micro Mill is designed to mill small bone pieces, for example from the mandible, into smaller bone chips. Due to its compact design, the Noviomagus Micro Mill can be hand held and requires little effort to mill bone into smaller bone chips.

Bone chips are collected in the Micro Chip Cup which can be unscrewed after the milling process to retrieve the milled bone. When finished, the Noviomagus Micro Mill can be fully disassembled with the Micro Wrench and all the separate parts can be autoclaved.



There are two types of milling drums for the Noviomagus Micro Mill. The *Extra Extra Fine* Micro Milling Drum produces chips between 1 and 2 mm. The *Extra Extra Extra Fine* Micro Milling Drum produces bone chips smaller than 1 mm.

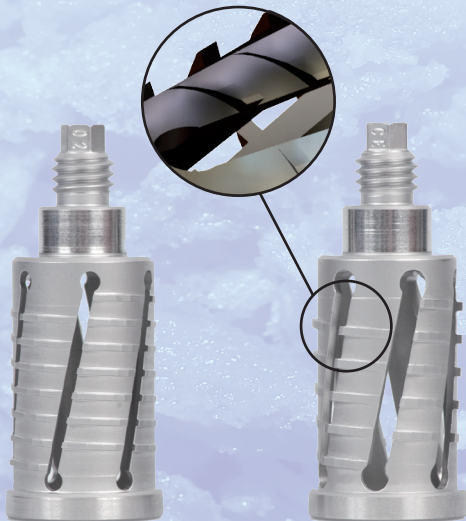

Two types of milling drums for various chip sizes. Each milling drum has spirally positioned teeth.

Advanced benefits

- Compact and intuitive design
- Easy to assemble
- Two types of milling drums:
 - Extra Extra Fine produces 1 - 2 mm chips
 - Extra Extra Extra Fine produces chips < 1 mm
- Easy to clean
- Easy collection of bone chips in the Micro Chip Cup

Today's technology guarantees high quality chips


Ordering information

Noviomagus

Micro Mill

Noviomagus Micro Mill:


| | |
|---|-------------|
| Micro Housing | SBM-I011-38 |
| Micro Push Block | SBM-I011-39 |
| Micro Handle | SBM-I011-40 |
| Micro Chip Cup | SBM-I011-63 |
| Micro Wrench | SBM-I011-66 |
| Micro Milling Drum Extra Extra Fine | SBM-I011-65 |
| Micro Milling Drum Extra Extra Extra Fine | SBM-I011-42 |
| Sterilization Tray (NCM) | SBM-I011-60 |

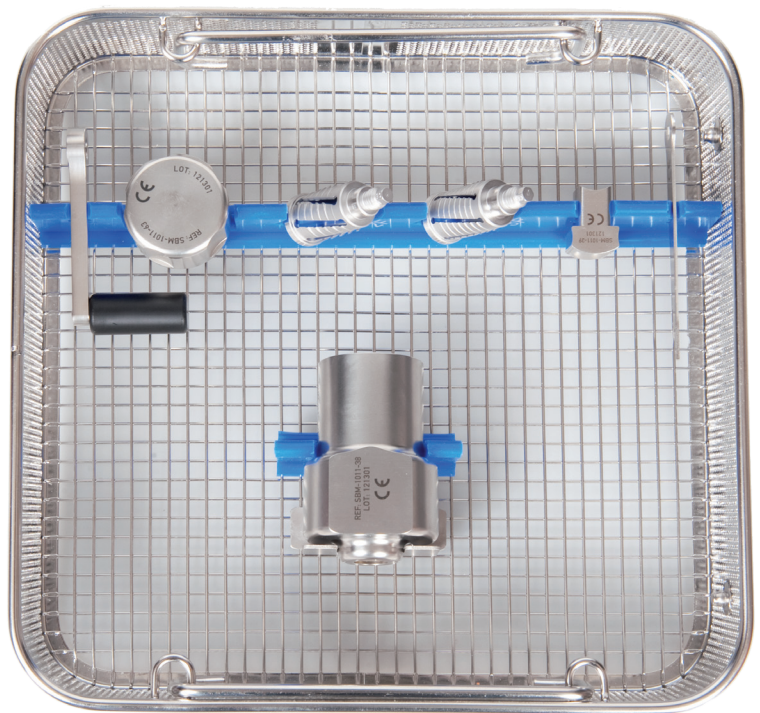

The Noviomagus Micro Mill is CE marked and FDA approved. US Patent No. 6,318,651. Other Noviomagus products manufactured by Spierings Orthopaedics B.V. are: Noviomagus Bone Harvesting System, Noviomagus Mini Mill, Noviomagus Bone Mill and Noviomagus Femoral Head Reamer Set.

No rights can be derived from the information in this brochure.

Copyright: Spierings Orthopaedics B.V.

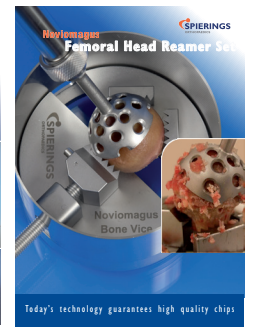

Distributed by:
HB-Medical B.V.
Tivolilaan 205
6824 BV Arnhem
The Netherlands
T: +31 -(0)26-3618804
F: +31 -(0)26-3882551
E: info@hb-medical.nl
I: www.hb-medical.nl