

Flexible Osteotome System

Provides an assortment of osteotome blades for various orthopedic surgery procedures

PRODUCT NO'S:

S0011-00	[Set with Quick-Coupling Handle and Case]
S0012-00	[Set with Locking Nut Handle and Case]
Sets Include / Available Individually:	
S1002	[Thin Osteotome Blade] 2.5" (6,3 cm) x 8 mm
S1003	[Thin Osteotome Blade] 2.5" (6,3 cm) x 10 mm
S1004	[Thin Osteotome Blade] 2.5" (6,3 cm) x 12 mm
S1005	[Thin Osteotome Blade] 2.5" (6,3 cm) x 20 mm
S1006	[Curved Thin Osteotome Blade] 2.5" (6,3 cm) x 12 mm
S1007	[Curved Thin Osteotome Blade] 5" (12,7 cm) x 20 mm
S1008	[Thin Osteotome Blade] 5" (12,7 cm) x 10 mm
S1009	[Thin Osteotome Blade] 5" (12,7 cm) x 8 mm
S1020	[Handle with Quick-Coupling End] 5" (12,7 cm)
or	
S1021	[Handle with Locking Nut] 5" (12,7 cm)
S1133	[Radial Osteotome] 5" (12,7 cm) x 10 mm
S1120	[Radial Osteotome] 5" (12,7 cm) x 12 mm
S1134	[Radial Osteotome] 5" (12,7 cm) x 14 mm
S1121	[Radial Osteotome] 5" (12,7 cm) x 16 mm
S1122	[Radial Osteotome] 5" (12,7 cm) x 20 mm
S2007	[Slap Hammer] 12" (30,5 cm)
9018	[Case]

Blade lengths reflect the actual working portion of the blade only. For overall length, add 1.5" (3,8 cm) to blade length listed above.



Handle with Quick-Coupling End #S1020



System Includes Choice of Handle Style


Handle with Locking Nut #S1021



2.5" (6,4 cm) x 8 mm Thin Osteotome Blade #S1002



2.5" (6,4 cm) x 10 mm Thin Osteotome Blade #S1003

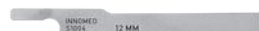

2.5" (6,4 cm) x 12 mm Thin Osteotome Blade #S1004



2.5" (6,4 cm) x 20 mm Thin Osteotome Blade #S1005



2.5" (6,4 cm) x 12 mm Curved Thin Osteotome Blade #S1006



5" (12,7 cm) x 20 mm Curved Thin Osteotome Blade #S1007



5" (12,7 cm) x 8 mm Thin Osteotome Blade #S1009



5" (12,7 cm) x 10 mm Thin Osteotome Blade #S1008



5" (12,7 cm) x 10 mm Radial Osteotome #S1133



5" (12,7 cm) x 12 mm Radial Osteotome #S1120



5" (12,7 cm) x 14 mm Radial Osteotome #S1134



5" (12,7 cm) x 16 mm Radial Osteotome #S1121



5" (12,7 cm) x 20 mm Radial Osteotome #S1122



Slap Hammer #S2007

- ▶ Sharp, flexible blades are well suited for loosening implants from cement or bony ingrowth fixation
- ▶ Various blade widths and profiles allow great flexibility to follow the implant contours
- ▶ Modular handle is made of high impact surgical stainless steel and has a quick-coupling positive locking mechanism for ease of use and quick blade changes
- ▶ Slap hammer threads into the handle and is designed to facilitate blade removal
- ▶ Optional Strike Plate can be attached to the Handle for direct striking with a mallet
- ▶ Optional Curved Chisel Blades are designed to help loosen the cement/prosthesis interval in TKA tibial tray and femoral component revisions. The curved design is useful in working around pegs & fins to get posterior cement access. Also helpful with removal of other implants, i.e shoulder, ankle, etc.

www.innomed.net



INNOMED

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com


1.800.548.2362

Flexible Osteotome System – Options

Optional Parts and Blades

PRODUCT NO'S:	
S1020-SP	[Strike Plate for Handle] Diameter 1.625" (4,1 cm)
Optional Osteotome Blades (Not Included In Complete Set):	
S1010	[Offset Osteotome Blade – 15 mm] 3.15" (8 cm) x 15 mm Offset Section Length: 1.1" (2,8 cm) Offset: .65" (1,6 cm)
S1123	[7.5" XL Osteotome Blade – 8 mm] 7.5" (19,1 cm) x 8 mm
S1135	[Radial Osteotome Medial Curve] 6.75" (17,1 cm) x 11 mm
S1136	[Radial Osteotome Lateral Curve] 6.75" (17,1 cm) x 11 mm
S1137	[Radial Osteotome Medial Curve] 5" (12,7 cm) x 11 mm
S1138	[Radial Osteotome Lateral Curve] 5" (12,7 cm) x 11 mm
Optional Chisel Blades (Not Included In Complete Set):	
S1233-L	[2" Flexible Left Curved Chisel] 2" (5,1 cm) x 8 mm
S1233-R	[2" Flexible Right Curved Chisel] 2" (5,1 cm) x 8 mm
S1222	[2.5" Chisel Blade – 8 mm] 2.5" (6,4 cm) x 8 mm
S1223	[2.5" Chisel Blade – 10 mm] 2.5" (6,4 cm) x 10 mm
S1224	[2.5" Chisel Blade – 12 mm] 2.5" (6,4 cm) x 12 mm
S1225	[2.5" Chisel Blade – 20 mm] 2.5" (6,4 cm) x 20 mm
S1229	[5" Chisel Blade – 8 mm] 5" (12,7 cm) x 8 mm
S1228	[5" Chisel Blade – 10 mm] 5" (12,7 cm) x 10 mm
S1231	[5" Chisel Blade – 12 mm] 5" (12,7 cm) x 12 mm
S1230	[5" Chisel Blade – 20 mm] 5" (12,7 cm) x 20 mm
S1227	[5.5" Long Chisel Blade – 8 mm] 5.5" (14 cm) x 8 mm
S1232	[7.5" XL Chisel Blade – 8 mm] 7.5" (19,1 cm) x 8 mm
S1234	[8.5" XL Chisel Blade – 8 mm] 8.5" (21,6 cm) x 8 mm
S1235	[9.5" XL Chisel Blade – 8 mm] 9.5" (23,1 cm) x 8 mm
S1236	[10.5" XL Chisel Blade – 8 mm] 10.5" (26,7 cm) x 8 mm
S1237	[11.5" XL Chisel Blade – 8 mm] 11.5" (29,2 cm) x 8 mm
S1238	[12.5" XL Chisel Blade – 8 mm] 12.5" (31,8 cm) x 8 mm

Blade lengths reflect the actual working portion of the blade only.
For overall length, add 1.5" (3,8 cm) to blade length listed above.

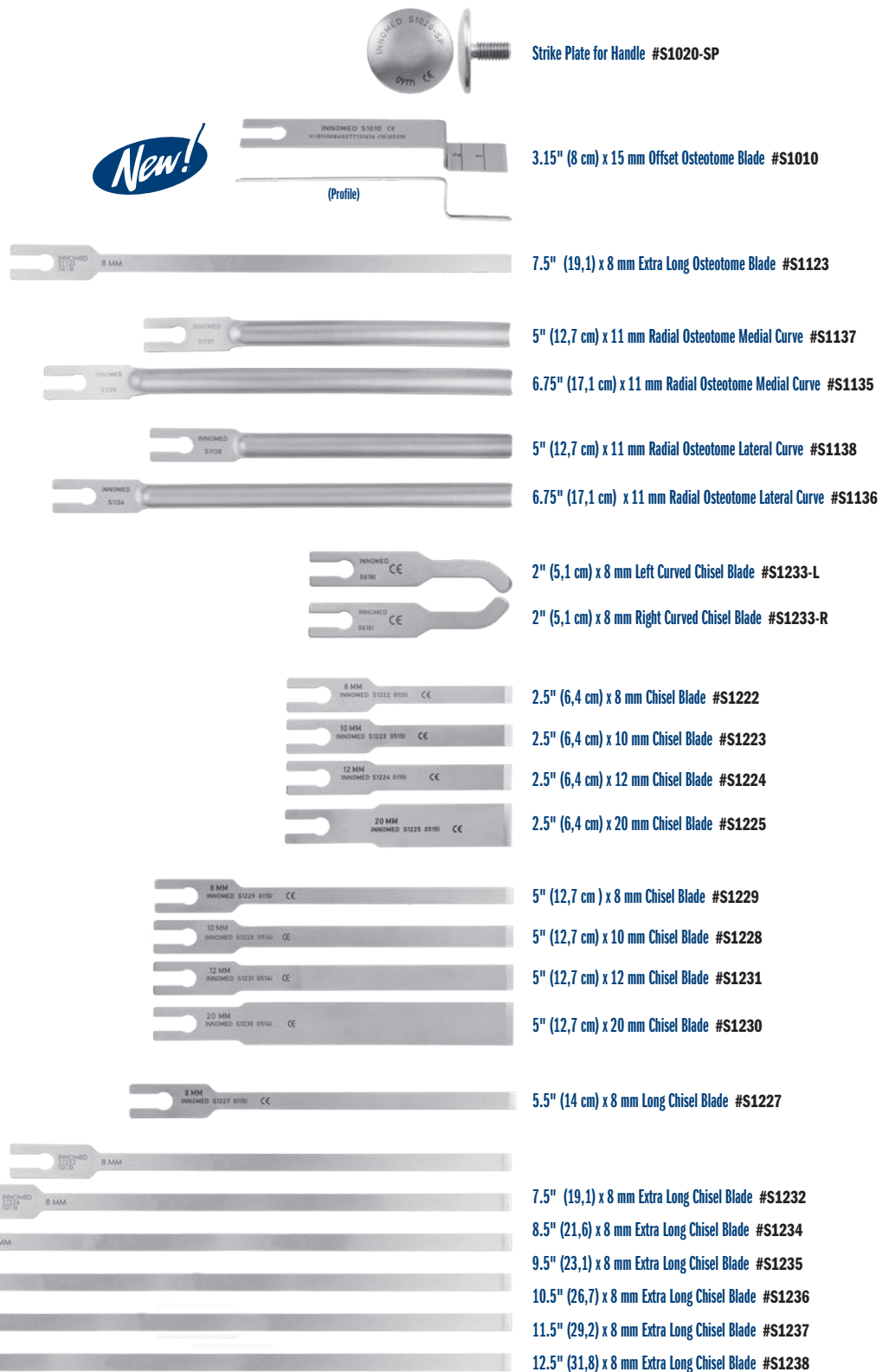

Extra Long Chisel Blades are designed for removal of well-fixed long bone intramedullary hardware

Medial and Lateral Curve Radial Blades designed by Henry Boucher, MD
Curved Chisel Blades designed by William McMaster, MD

www.innomed.net



INNOMED

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362